


Basics on Biosignal Measurement with g.MOBIIlab⁺ V2.14.00


Unipolar / Bipolar Biosignal Derivations

Connecting External Sensors to g.MOBIIlab


Unipolar EEG recording with g.MOBilab

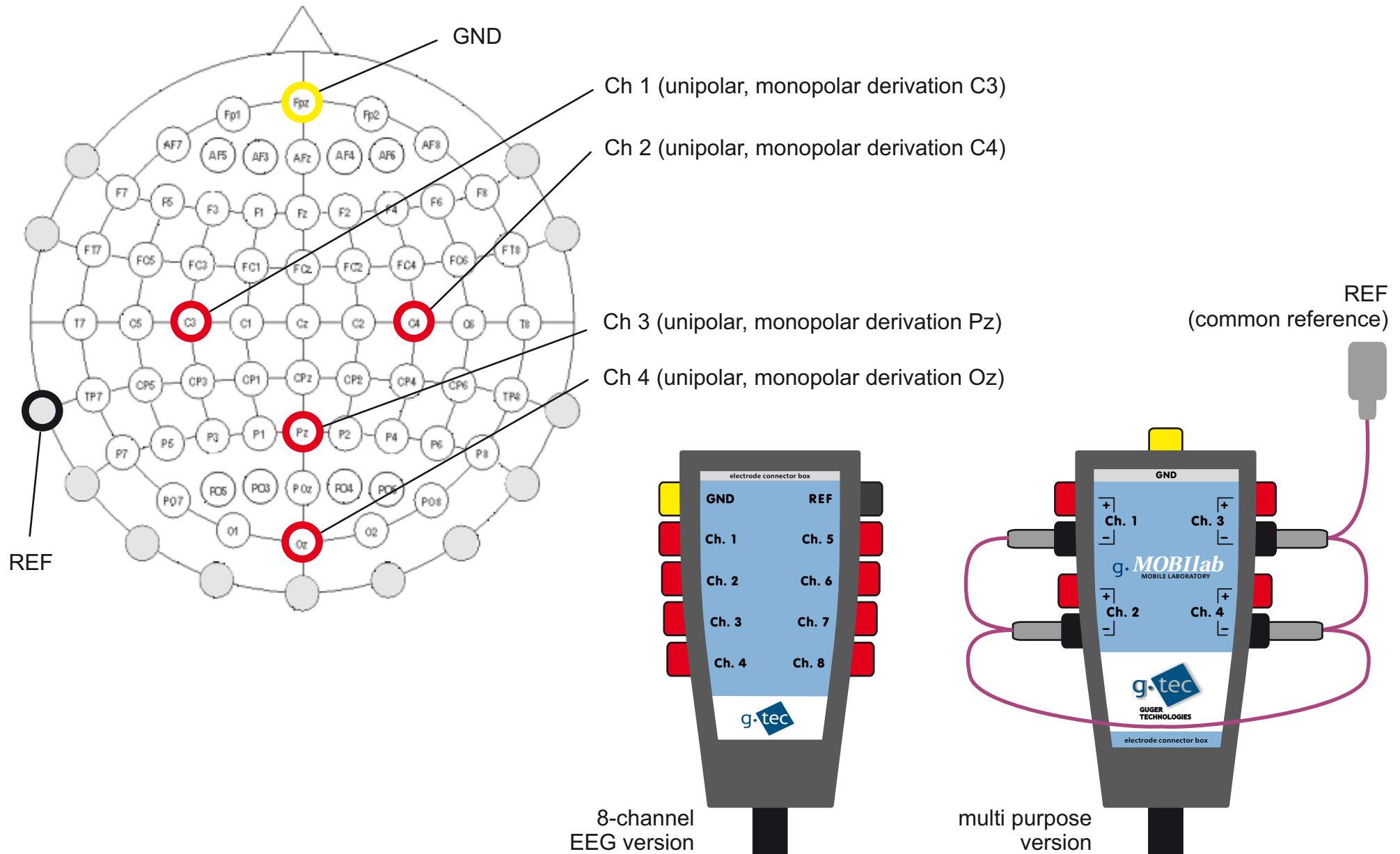

Bipolar EEG/EOG recording with g.MOBilab

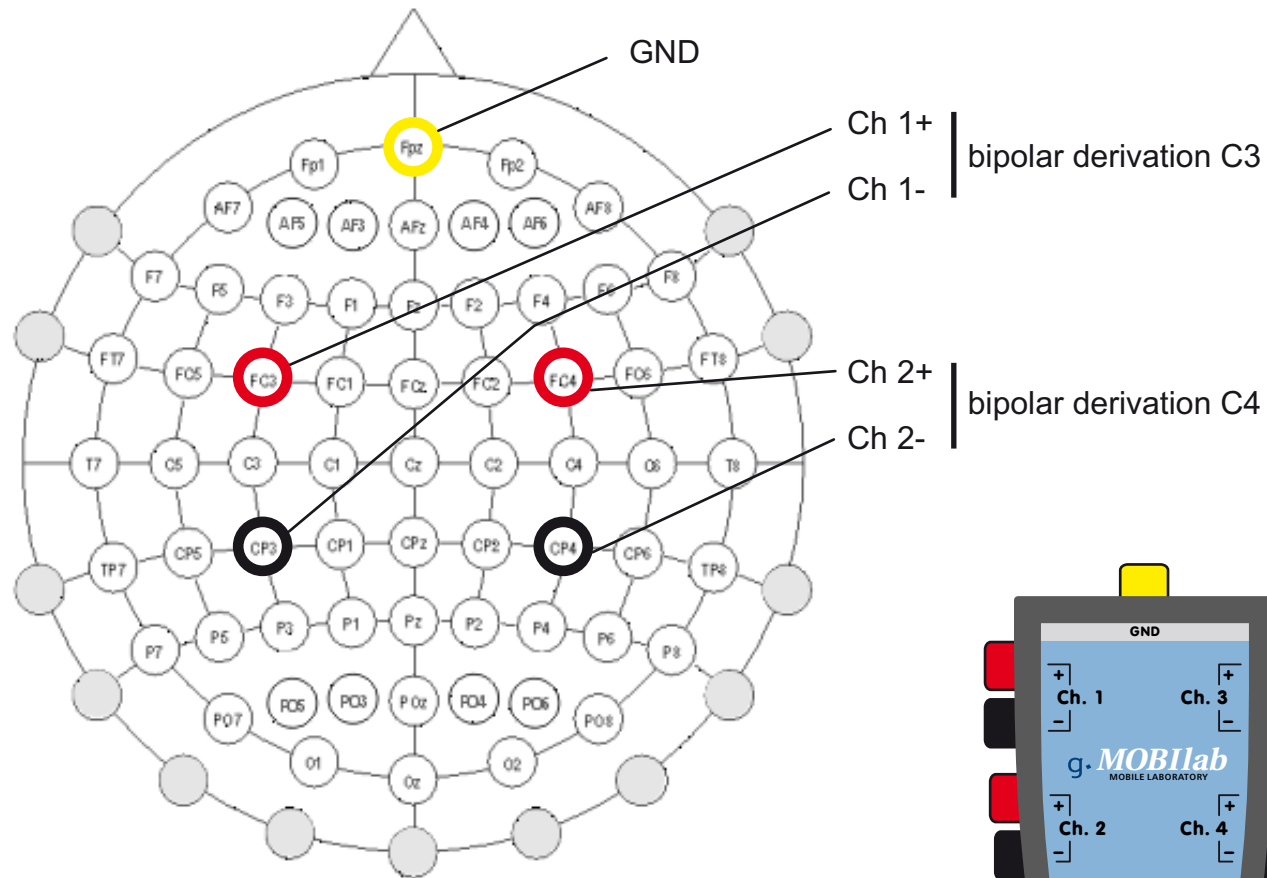


GND
(only required
if not connected
to Fpz)

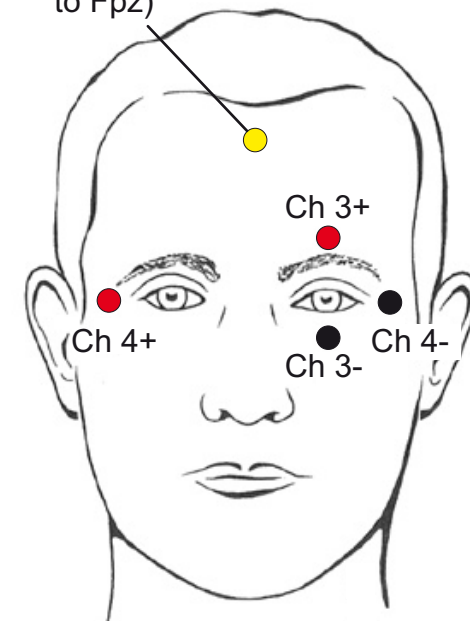


In this example:
Ch 3: VEOG
Ch 4: HEOG



multi purpose
version

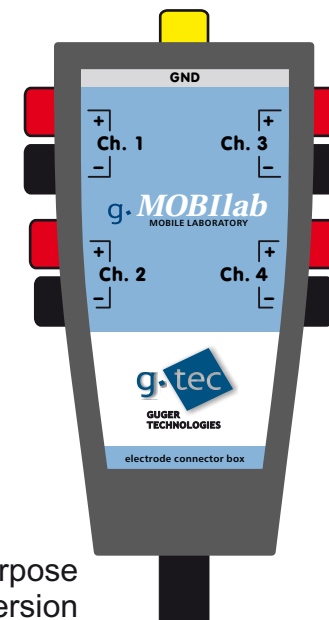

Bipolar ECG/EMG recording with g.MOBIIlab

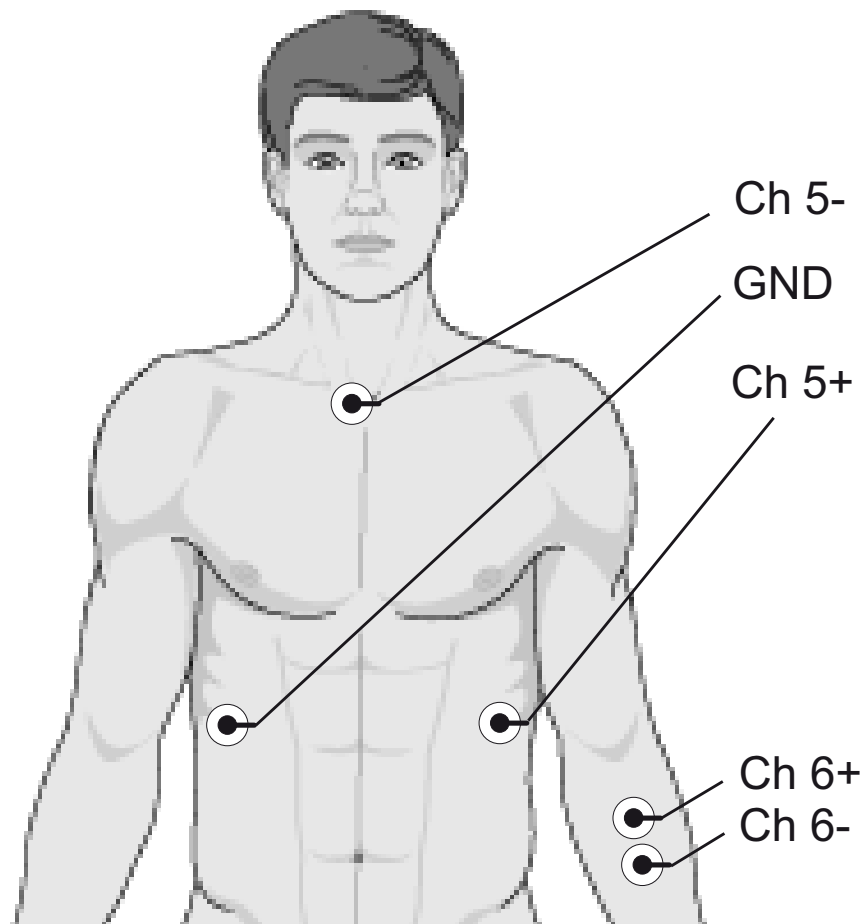


example: Ch 5: ECG, Ch 6: EMG

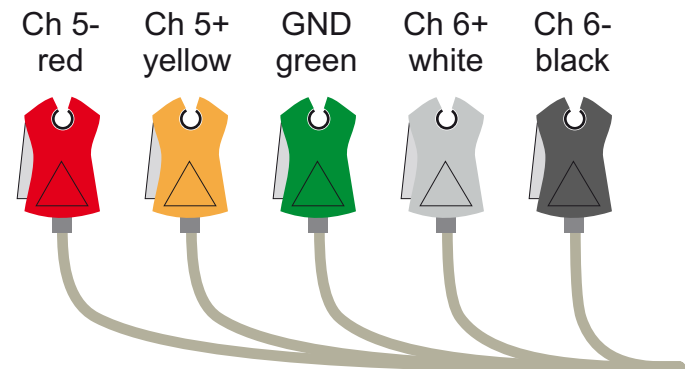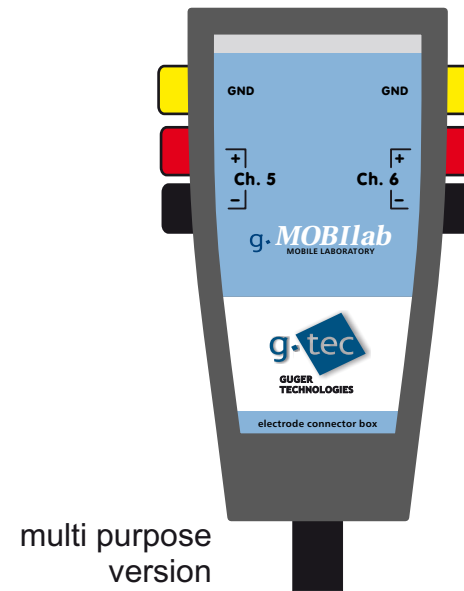

5-lead ECG/EMG patient cable

Connecting external sensors to g.MOBllab

Sensors with a 2-lead sensor cable
(e.g. snoring microphone, respiration sensor, ...)

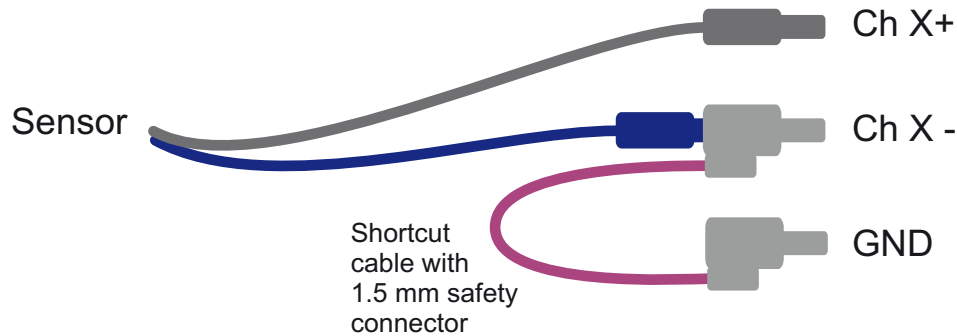

Sensor type:

Recommended input channel:

g.PULSEsensor
(low output)

Ch 3 or Ch 4
(multi purpose version)

g.PULSEsensor
(high output)

Ain 1 (Ch7) or Ain 2 (Ch8)
(multi purpose version)

g.RESPIRATIONsensor

Ch 3 or Ch 4
(multi purpose version)

g.SNORINGSensor

Ch 1, 2, 5 or 6
(multi purpose version)

g.GSRsensor

Ain 1 (Ch7) or Ain 2 (Ch8)
(multi purpose version)

g.TEMPsensor

Ain 1 (Ch7) or Ain 2 (Ch8)
(multi purpose version)

Sensors with a 3-lead sensor cable
(e.g. g.PULSEsensor)



Sensors with DC-high level output
(e.g. g.GSRsensor, g.TEMPsensor)



Anti-Static Kit

Avoid or reduce artifacts in biosignal recordings resulting from electro-static charges in a sub-optimal lab environment and protect sensible electronics.

