

How to install and setup the Bluetooth connection of g.MOBllab+ for paired Bluetooth devices, Version 3.14.01


The purpose of this tutorial is to instruct the user to set up and configure the paired Bluetooth connection to g.MOBllab+ based on the Ezurio Bluetooth USB dongle. This applies to Microsoft Windows 7 Professional English Bluetooth driver stacks. The device pairing has only to be done when the device is the first time installed or anytime you uninstalled g.MOBllab+.

1. Insert the Bluetooth dongle into a free USB slot. Wait until the Microsoft Bluetooth driver is installed automatically and the Bluetooth icon appears in the Notification area in the bottom right corner of your screen.




Double-click this icon.


2. Make sure that g.MOBllab+ is switched on and the green LED is blinking.
3. In the **Bluetooth Devices** dialog click on the **Add a device....**

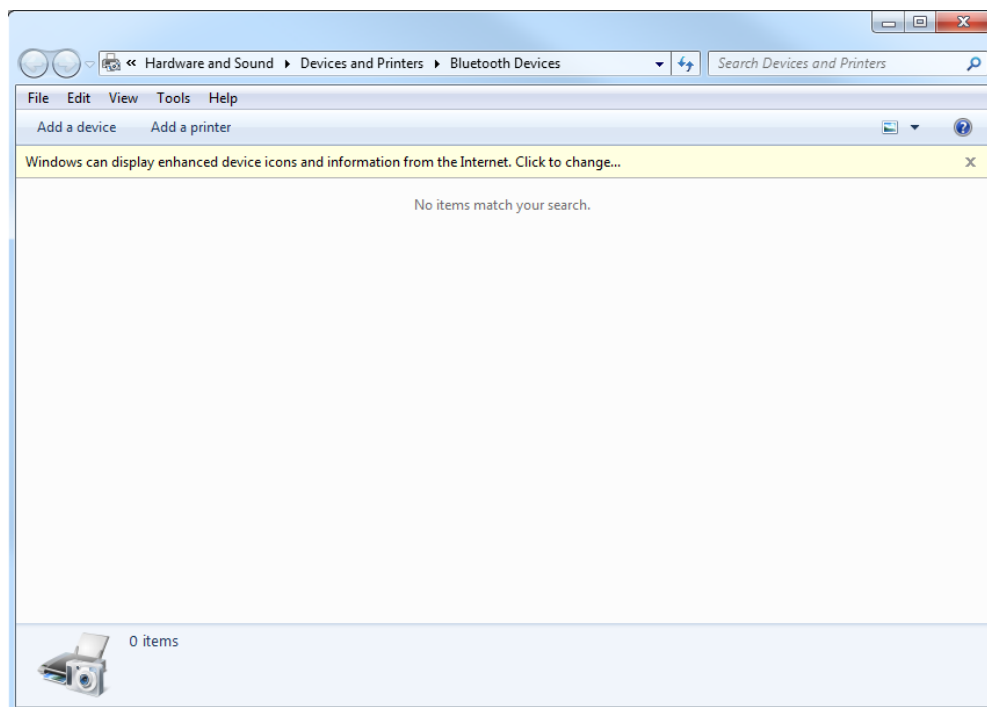

4. Select the entry with the g.MOBllab+ serial number (e.g. MP-2009_01_14) and click **Next**.

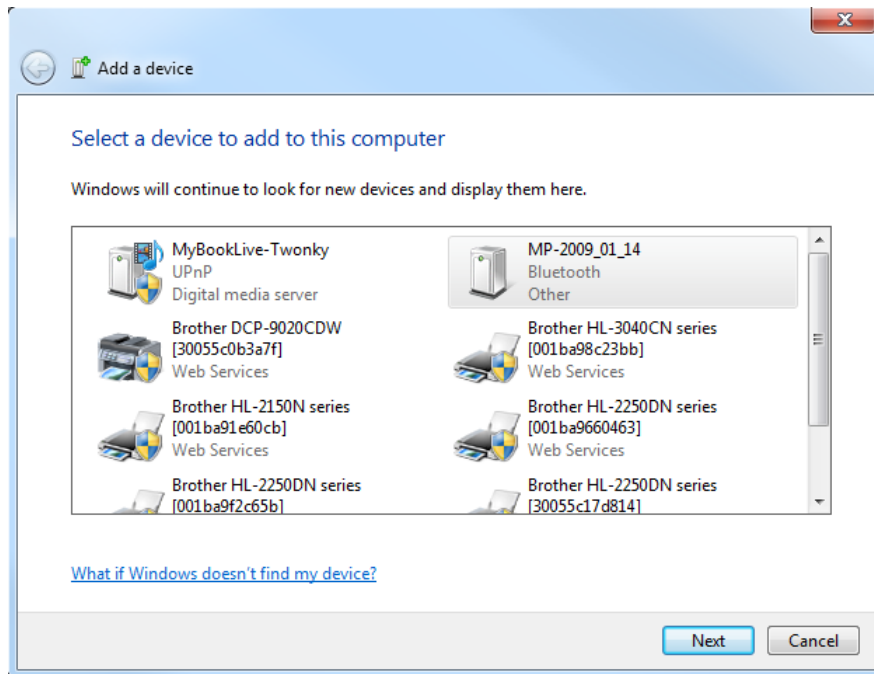


5. A passkey is needed to establish a connection to g.MOBllab+, which can be found in the battery compartment of g.MOBllab+ and is normally 1234. Enter 1234 as passkey and click **Next**.

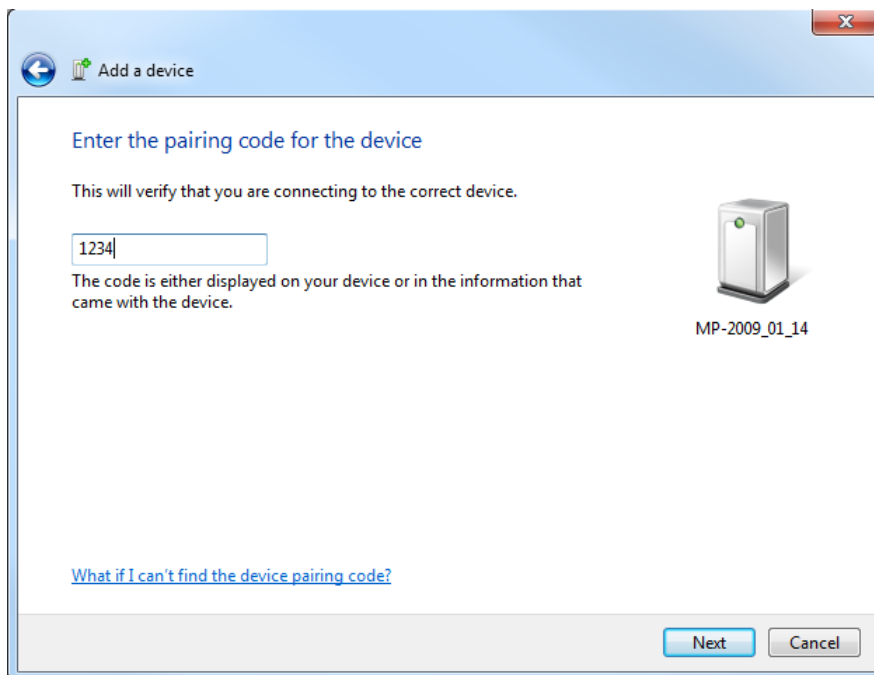

6. When installation is complete, the dialog below is shown. Click **Close** to exit the wizard.



7. To get the COM port where g.MOBllab+ is connected, right-click **MP-2009_01_14** and select **Properties**.

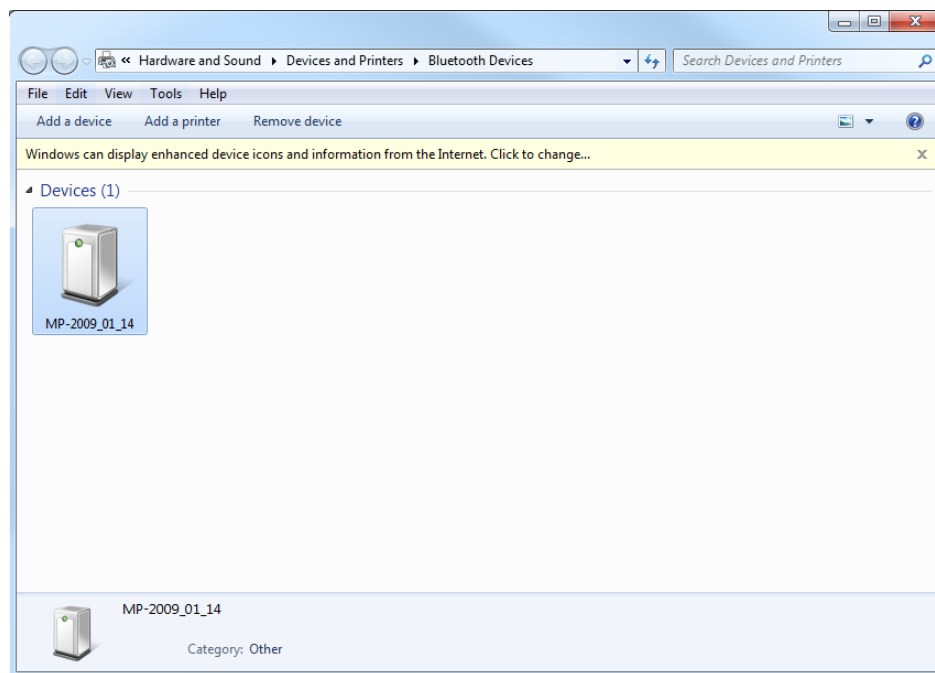

8. Select the **Services** tab, there is the serial port of g.MOBllab+, **Serial port (SPP) 'AT Serial'**, in this case **COM5**.

